UMINATE ME

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size : 16.5cm x 12.5cm x 4.5cm

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The logo consists of a 4x4 grid of squares. The top-left square is white, the top-right square is black, and the remaining 15 squares are white.

UMINATE ME

User Manual

Serial Number: UM-V2

Production : August 2024

Introduction

The **Dahan Tab** is designed to help accelerate tooth movement during orthodontic treatment and is used intraorally.


This device emits light with a wavelength close to infrared, produced by a specialized LED.

This wavelength can trigger a suitable cellular response in the tissue surrounding the tooth following orthodontic force application, and it can increase the rate of tooth movement by up to 36%.

The device consists of two parts: an intraoral section and a control unit. The control unit contains a 3.7-volt lithium battery, which is charged using the device's accompanying charger.

2

Intra Oral Part

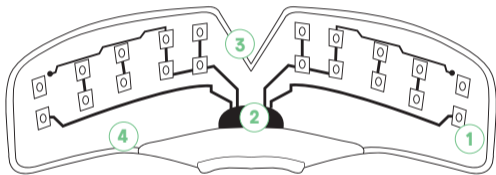1

LEDs

4

Body

2

USB Port

3

Anterior

3

Control Part


1_ Notification

2_ On/Off

3_ USB Port

4_ Charging Port

5_ Mode Button

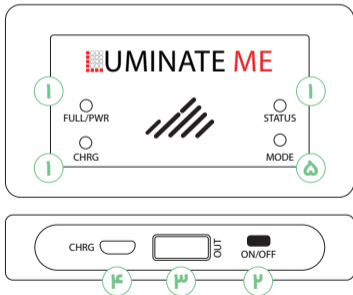Check List

Manual Guide

Charger

5 Device Preparation and Setup

Turn on the device using the ON/OFF button until the green LED on the top panel lights up. At this point, the device is ready for the user to set the operating mode.

- Mode One: In this mode, the indicator light blinks once and the device emits a sound. Then, all LEDs turn on and operate continuously.
- Mode Two: In this mode, the indicator light blinks twice and the device emits a sound. Then, all LEDs turn on and operate intermittently.
- Mode Three: In this mode, the indicator light blinks three times and the device emits a sound, but none of the LEDs turn on. This mode is used for research purposes.

6

Device Usage

The user has 5 seconds to change the mode each time the mode button is pressed. If the user does not change the mode within this time, the device will emit a long beep sound and enter the selected mode, remaining active in this mode for 300 seconds.

You can now place the oral module inside your mouth, gently position it with your teeth, and keep it in the mouth for the default time or the duration recommended by your dentist.

After the set time has elapsed, you can remove the device. The device will emit a beep sound for 5 seconds to indicate that it has completed its task, and the green LED will remain steady. The device will then be ready for reprogramming and setting.

7

Health Tips

After each use, detach the oral section of the device from the USB connection and clean it with a toothbrush, toothpaste, and cold water.

When cleaning the device, make sure that water does not reach the USB section.

Please note that during the initial sessions of using the device, increased saliva production is normal. However, if this issue persists beyond the third session, please inform your treating dentist or the product support team.

-  The device's LED emission is in the infrared spectrum and is not visible under natural light. Avoid prolonged exposure to it.
-  Since infrared LED emission may not be easily detectable in bright environments, if you need to check its operation, do so in a relatively dark space or use a mobile camera.
-  Avoid connecting the device to other electronic devices with USB inputs.
-  If the groove at the front of the device is shaped like a "V," it can be used for the upper jaw. If the groove is shaped like an "Λ," the device is suitable for use with the lower jaw.

9

Reuse

Before reuse, ensure that the device, especially the USB section, is completely dry. Then reconnect it to the external oral section and proceed with using it.

It is necessary to fully charge the device using the accompanying charger after every two uses.

Since this device uses nickel-cadmium batteries, disposing of it in regular waste can cause environmental pollution. Therefore, after completing the treatment, please return the device to the company's address for recycling and take advantage of the offered incentives.

Contact Us

Now that you have finished reading the manual, you can begin using the device.

Please follow "**Dahan Tab**" on social media, and feel free to share any feedback or suggestions with us.

Feedback

lumosmed.co@gmail.com

Website

www.lumosmed.co

Support

☎ 0113842667733

Fara Noor Salamat

Hamedan, Khajeh Rashid Crossroad,
Health Technology and Innovation Center,
Hamedan University of Medical Sciences

081 - 3826 67 33
www.lumosmed.co

size : 16.5cm x 12.5cm x 4.5cm

Made in IRAN